


Ludger-Velocity™: A fast system for glycan labelling

Your Glycoprotein

Step 1:
Release glycans

PNGase F enzyme
E-PNG-01 or E-PNG-05

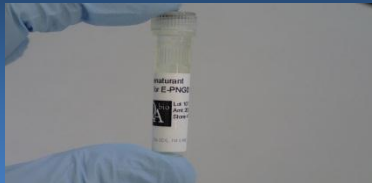

3 hours *

Dry down samples

2 hours

Step 2: Label

LT-KAA-VP24 kit
or
LT-KAB-VP24 kit



1 hour

Step 3: Clean

LC-T1-A6 cartridges
used with vacuum
manifold



1 hour

Sample is ready for analysis within a day

*It is advisable to perform a time course from 1-24 hours for PNGase F release